



# CD7

## EMERGENCY Evacuation Chair

See video demonstration  
[www.baronmead.com](http://www.baronmead.com)

**Baronmead**  
INTERNATIONAL

Established 1982

**Tel: +44 (0)1243 586692 Fax: +44 (0)1243 586312**

Flansham Business Centre, Hoe Lane, Flansham, Bognor Regis, West Sussex, England. PO22 8NJ

Web: [www.baronmead.com](http://www.baronmead.com) E-mail: [sales@baronmead.com](mailto:sales@baronmead.com)




- ◆ Passengers weight moves the unit downstairs, while the governor brake controls speed
- ◆ Failsafe brake brings the unit to a complete stop automatically
- ◆ Adjustable safety straps
- ◆ Stable and self-supporting
- ◆ Unique seat design allows easy transfer from wheelchair
- ◆ Durable track grips stairs securely
- ◆ Easily stores in a secure area when not in use
- ◆ Requires minimal maintenance
- ◆ Full aftersales service, parts & training



## Technical Specifications

Capacity	130kg (287lb/20.4st)
Dimensions (h x w x l) deployed	810mm x 426mm x 1310mm (32" x 16.8" x 51.5")
Dimensions (h x w x l) stored	1100mm x 426mm x 270mm (43.3" x 16.8" x 10.5")
Weight	20kg (44lb)
Max. stair angle	40°
Speed	1.1m/sec. (varies with weight/angle)
Speed control	Via hydraulic retarder
Brake	Operator controlled drum brake

Passenger/operator weight equivalence + / - 50% (recommendation only)

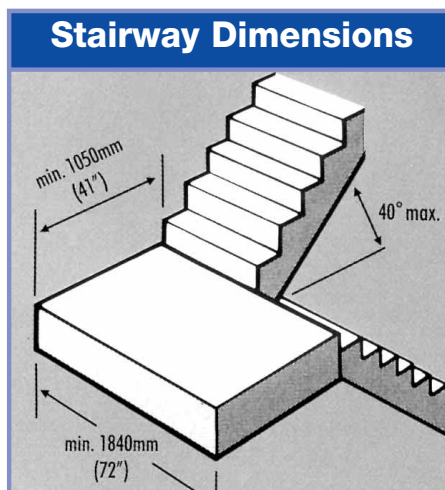

**Baronmead**  
INTERNATIONAL

Established 1982

Flansham Business Centre, Hoe Lane, Flansham,  
Bognor Regis, West Sussex, England. PO22 8NJ

**Tel: +44 (0)1243 586692**

Fax: +44 (0)1243 586312

Web: [www.baronmead.com](http://www.baronmead.com)

E-mail: [sales@baronmead.com](mailto:sales@baronmead.com)